

Introducing the



HOVERMATT®

Single-Patient Use

Air Transfer System

The #1 choice of hospitals for lateral patient transfers is now available in a disposable option. Just like the reusable HoverMatt transfer mattress, the Single-Patient Use HoverMatt transfer mattress is ideal for use in any patient care department. With the HoverMatt transfer mattress product line...

- Eliminate injuries related to lateral transfers and repositioning
- Improve patient safety and comfort during transfers
- No weight limit






HOVERMATT®

Single-Patient Use

HM34SPU*

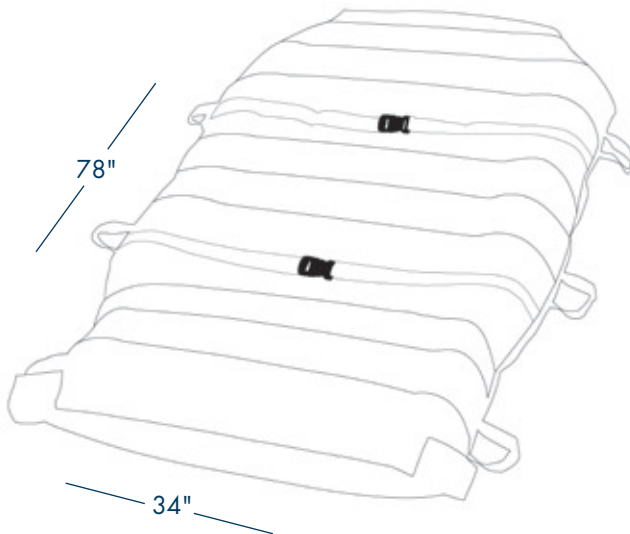

*Compatible with current HoverMatt Air Supply

Patent Pending

HM34SPU DIMENSIONS (INFLATED)

34" W x 78" H

HM34SPU PER PACKAGE

Box of 10

HM34SPU MATERIAL SPECS

Sewn Construction

Top: Non-Woven Polypropylene Fiber

Bottom: Nylon Fabric

Latex Free

Radiolucent/MRI compatible

Water Soluble Perimeter thread (120°) – Thread in sewn handles is NOT Water Soluble

The Single-Patient Use HoverMatt transfer mattress significantly reduces the possibility of cross contamination. The mattress is constructed with water soluble thread, which ensures that a contaminated product is not re-used – eliminating the need for outside laundering services. The mattress' white top will stain, therefore alerting the staff to potential infection concerns.

In the event the Single-Patient Use HoverMatt transfer mattress is laundered, it will no longer inflate when connected to the air supply; the "DO NOT WASH" label at the foot-end of the product will dissolve and state, "DISCARD IMMEDIATELY."



513 S. Clewell Street
Bethlehem, PA 18015

Phone 800.471.2776

Fax 610.694.9601

www.HoverMatt.com